

# CLASSIC KTP LASER

The STAR of Photocoagulation Systems

- > KTP Laser Module with two laser ports
- Advanced pneumatic micromanipulator enables SAFE and ACCURATE positioning
- Intuitive, single touch control dial selects treatment parameters
- ➤ Laser delivery and slit lamp illumination are ALWAYS coaxial



## PCL 5 ZL

#### SLIT LAMP & LASER

The laser optic of the PLC5 is unique because it enables coagulation in the periphery, which is a dramatic difference from other slit lamps where the illumination comes from below.

The laser optic is integrated into the slit illumination to assure coaxial viewing with laser delivery. The compact laser console can be positioned either to the right or the left of the slit lamp. Since the connecting cables are concealed, the surgeon can operate the laser without any restrictions.

The small laser package fits easily into any office setting. Portability can be enhanced by selecting a slit lamp table with wheels (optional).







### PCL 5 SHL


#### SLIT LAMP & LASER

The SUPER STAR of Photocoagulation Systems! Precision optics combined with a fully coaxial laser coupled to the slit lamp, sets the new standard in today's market. Unique features such as the neutral colour design filter and the noiseless micromanipulator combined with the hand trigger enables the

ultimate in performance of laser photocoagulation. With the laser console on the table, all functions are easily accessible by your finger tips from the base of the slit lamp.

It's what is known as **PERFECTION**.

PCL 5 SHL

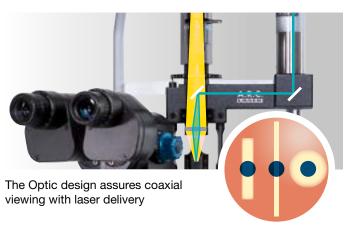



Single hand controls to trigger the laser pulse and adjust the aiming beam intensity. Suitable for left and right hand operation.





The Classic is not only a photocoagulator, but is a device that can be configured to best satisfy your needs.

Protection filters automatically switch in or out, which allows a perfect view when the laser is not is use.

With up to 1.5 Watts of power available to the cornea, it is possible to consider the Classic for multiple applications. The laser is compact, light weight and portable. The large display screen is highly visible from every angle.

Adaptors to a selection of slit lamps and operating microscopes are available.





An optional battery pack allows "cable free" use inside the OR – maximizing convenience.



Easy transport and stow away with the tailored transportation case.

Automatic or permanent Neutral Color filter design.



Choose from a wide range of contact lenses.





# CLASSIC & PCL5 ZL + PCL5 SHL

#### **Quality Made in Germany**

#### CLASSIC IS COMPACT AND EXTREMELY VERSATILE - ONE DEVICE, WITH MANY OPTIONS

- > Slit Lamp Adaptors
- > Microscope usage with a selection of Endo Probes
- > Laser Indirect Ophthalmoscope (LIO)
- > and more....

Our commitment is to supply our customers with the best quality products with the latest innovations and broadest range of accessories – available from a single source.

A.R.C. Laser in Germany is a privately held company to better serve you.

#### LIO



#### **ENDO PROBES**

A.R.C. Laser offers a selection of single use probes at an affordable price.





#### SPECIFICATIONS SLIT LAMPS

#### PCL5 SHL

Slit width	0 to 8 mm continuous
Slit length	1 to 8 mm continuous
Illumination field diameter	0.2, 1, 2, 3, 4, 6, 8 mm
Slit rotation	± 90° (0-180°)
Illumination unit swing	horizontal ± 90° vertical 0-20°
Filters	blue, redfree (green), grey (10%) UV and IR filter are built-in
Illumination	Halogen 6 V, 20 W 0 lux to 600,000 lux (Halogen)
	optional: LED 14 V, 5 W
Eye protection filter	permanent
Laser optic	parfocal 50 - 600 µm
Power requ.	110 V - 240 V 50/60 Hz
Weight	12.7 kg

#### PCL5 ZL

Slit width	0 to 16 mm continuous
Slit length	0.5 to 16 mm continuous
Illumination field diameter	0.2, 1, 2, 3, 4, 6, 8 mm
Slit rotation	± 90° (0-180°)
Illumination unit swing	horizontal ± 90°
Filters	blue, redfree (green), grey (10%) UV and IR filter are pre installed
Illumination	Halogen 6 V, 20 W 0 lux to 600,000 lux (Halogen)
	optional: LED 14 V, 5 W
Eye protection filter	permanent
Laser optic	parfocal 50 - 600 µm
Power requ.	110 V - 240 V 50/60 Hz
Weight	12.7 kg

Stereo microscope	
Stereo angle full	10°
Oculars	± 7 diopters
Magnification	6.3x / 10x / 16x / 25x / 40x
Pupillar distance	52 to 78 mm

#### **CLASSIC GREEN LASER**

Laser	KTP - 532 nm
Output power	1.5 Watt max.
Pulse length	10 ms to cw
Pulse interval	100 ms to 1,000 ms
Frequency	1 to 9 Hz
Aiming beam	635 nm, 1 mW max.
Adapters	Different slit lamps
Foot print	250 x 190 x 52 mm
Weight	3.3 kg - Complete
Power supply	Rechargeable Batteries or 90 to 230 VAC
Fiber Attachments	2 Fiber ports
Accessories	Indirect Ophthalmoscope / Endo Probes
* at fiber end The device and the delivery system is protected by one or more patents. Technical data may change without notice.	AVOID EXPOSURE TO BEAM CLASS a BLASER PROJUCT BUCK AND BEAM CLASS a BLASER PROJUCT BUCK BASER BASER PROJUCT BUCK BASER PROJUCT

Your local distributor:



A.R.C. Laser GmbH Bessemerstraße 14 D-90411 Nürnberg Germany t +49 (0) 911 217 79 -0 +49 (0) 911 217 79 99 info@arclaser.com www.arclaser.com